

Alert-it Care Systems


P137& P138 Care Call Pager System

User Handbook



One of a range of Alert-it Care
Alarms available from:

After 30 sec inactivity the display returns to STANDBY with backlight off.

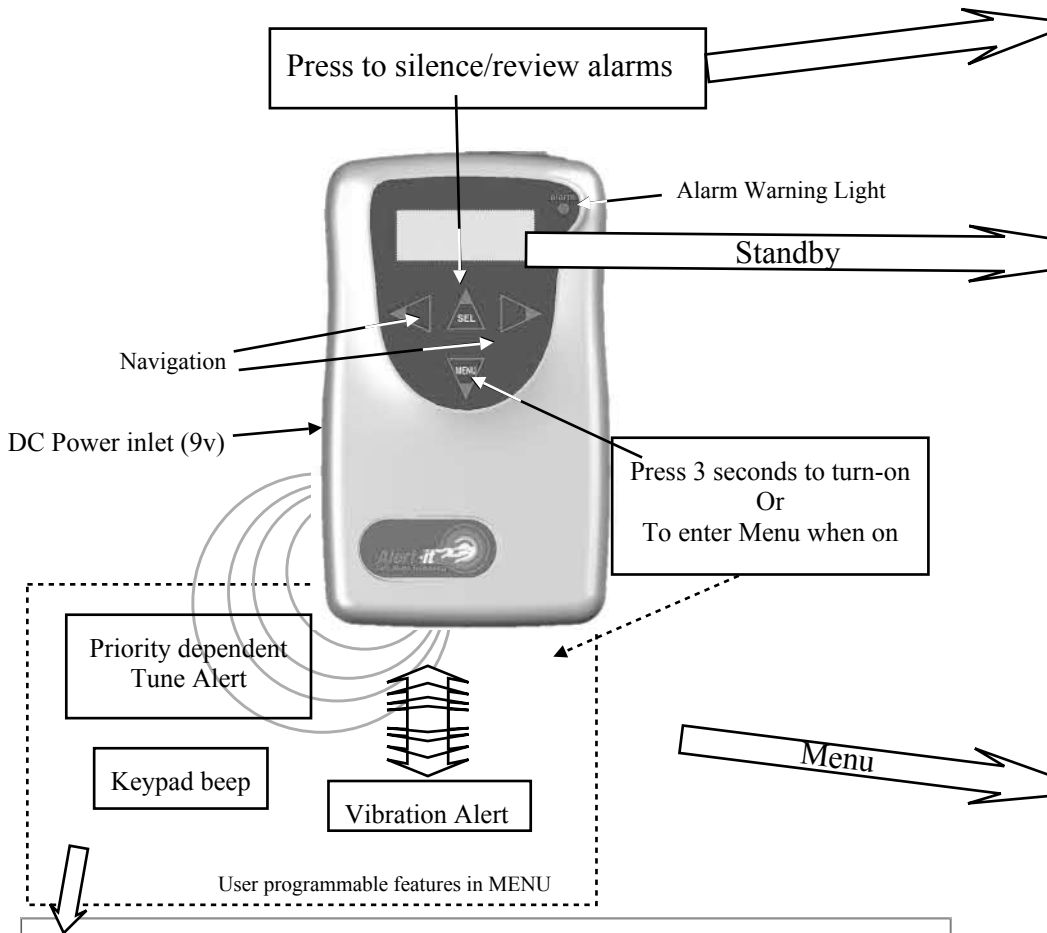

Silencing RF Fail Alarms.


The RF Fail alarm is used to warn of possible communication failure that would leave the client vulnerable. However it is also a natural consequence of turning off the client equipment, which may be required to prevent false alarms when the client is absent.

In this case there is the option to silence the tune that accompanies this alarm. The display will still warn that the equipment is off, but the frequent audible alarm will no longer cause nuisance. Any true distress alarm will over-ride the display and alert the carer as required.



The Alert-it system has been designed with due regard to reliability and integrity. While it offers a highly vigilant monitoring method, it is always possible that a distress condition can go undetected for a variety of reasons (including malfunction) and in life threatening situations it is advisable to use the Alert-it system in conjunction with additional monitoring techniques (e.g. video). Neither the manufacturer nor its agent can accept legal responsibility to provide a system that is infallible.



◀ Navigate all active alarms ▶

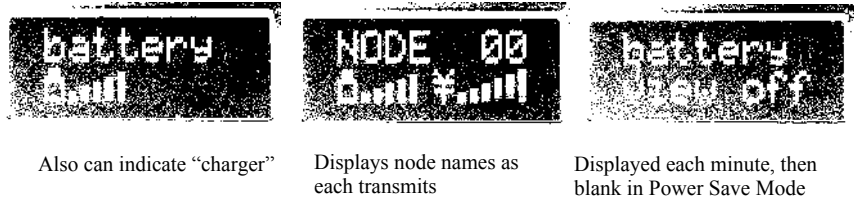

Optional stand-by displays

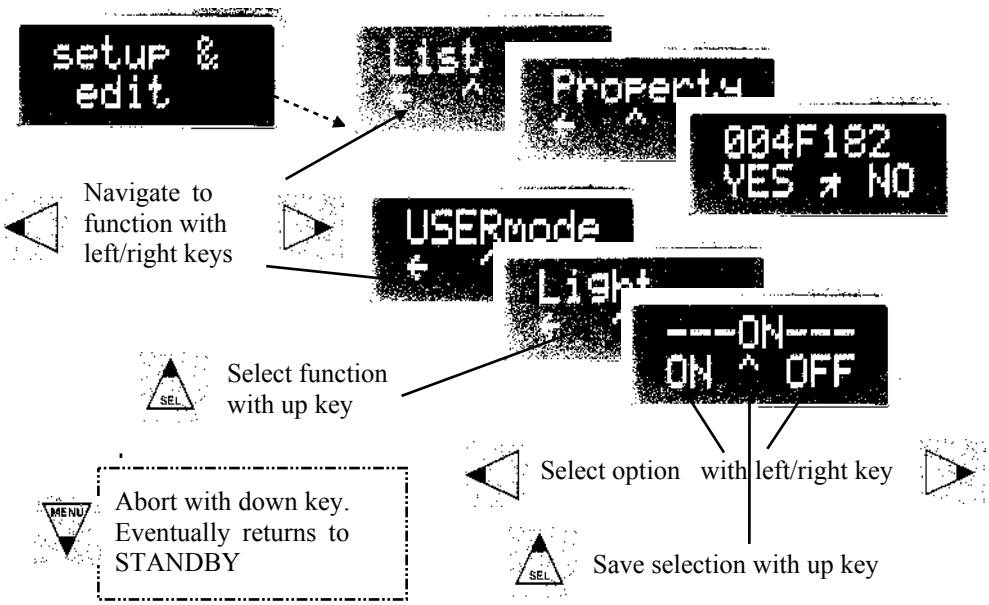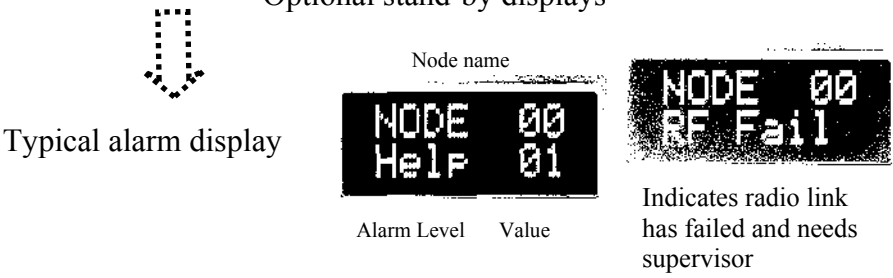


Extended User Functionality

The more advanced features of the Pager are barred from the normal user by a 4 digit password.

The supervisor can extend the user functionality to cover some of these advanced features such as:

- Turn the pager off
- Reset all alarms
- Modify the Name associated with a node

For details of these features please refer to the Installation Handbook UH1068A

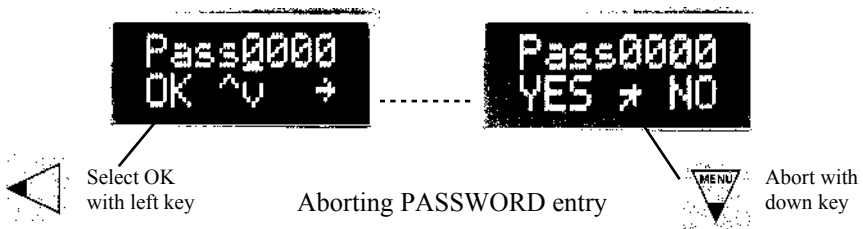

Battery Maintenance

The Pager is powered by two rechargeable AA cells. These are charged in situ if:

- A 9v adapter (1.3mm pin) is connected to the DC inlet socket
- Or the unit is in the P146 USB Interface with a 12v adapter (2.1mm pin) connected.

With power applied the display will indicate “charger” in standby, instead of “battery” with two different icons for fast and standby charge.

On fast charge the pager may feel warm to touch, which is normal.

If, after prolonged use, the batteries need replacement, then open the rear compartment by sliding the cover down, and replace with NiMH batteries only.