

This system is certified to the following European Standards

93/42/EEC¹	Class 1 Medical Device
73/23/EEC	Low Voltage Safety
BS EN 61010-1	Radio Interference Immunity
89/336/EEC	
BS EN 50081-1:1992	Permitted radio transmission Permitted Materials
BS EN 50082-1:1995	
BS EN 300-220-1	
2002/95/ECRoHS	

¹Alert-it Care Alarms are social aids designed and manufactured in accordance with 93/42/EEC as Class 1 Medical Devices. They are intended to improve the vigilance of carers to distressing side-effects of various diseases, such as Epilepsy and Dementia. They do not monitor vital physiological processes and should not be expected to diagnose any disease or predict the onset of any symptoms.

Support

For technical support please fax or EMail:
HELP: 0845 2179951
FAX : 0845 2179953
Support@itsdesigns.co.uk

Designed by:
ITs Designs Ltd
Leicester
LE9 9FE UK

...using technology to care for carers



None of the components, including batteries should be disposed of as Domestic Waste. For information on disposal contact ITs Designs Ltd.

Alert-it Care Systems

P156 Call Button

User Handbook



One of a range of Alert-it Care Alarms available from:


The P156 Call Button comes in 2 variances.

P156A has a 6 pin connector and is used with the P133G input extender to allow up to three switched sensors (eg a Bed Occupancy Mat and Call Button) to share the single Input C connector on the P154 Bedside Monitor. The P133G plugs into the 6way socket at the end of the P157 Extension lead, which in turn plugs into Input C.

The P156B has a larger 8 pin connector with a blue outer. This plugs directly into P154 Input A when no Bed Movement sensor is required, thus freeing Input C for the lone use of the Bed Occupancy sensor.

Pressing the orange area will generate an immediate alarm. The alarm is reset by pressing the RESET button on the P154.

On the P137/P138 Pagers the alarm will be presented as
URGENT 11

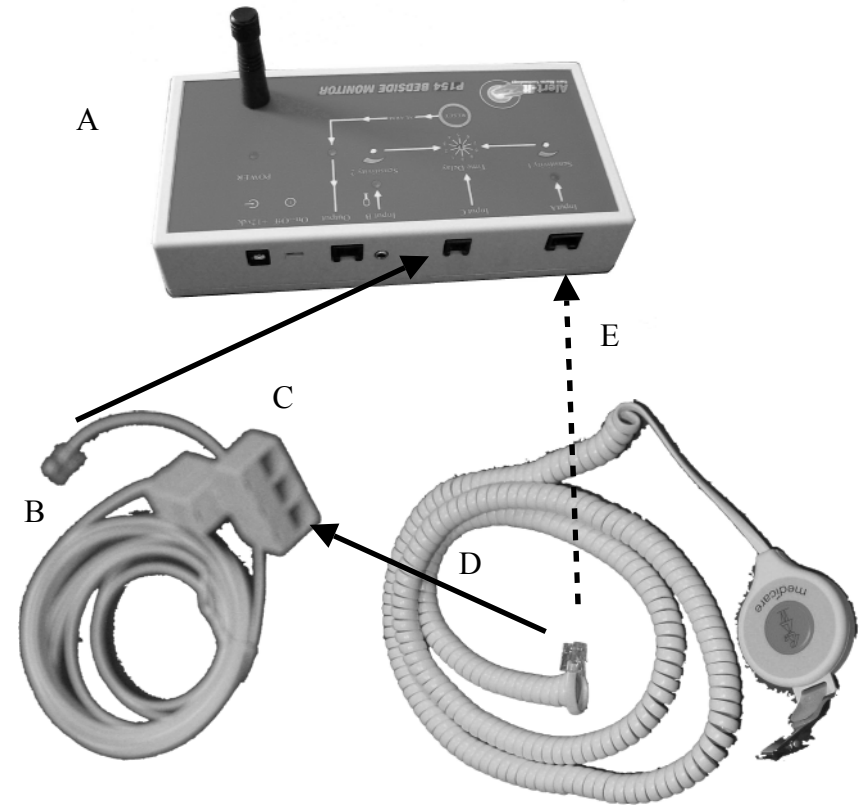

Part Description	Part No	Handbook
A Bedside Monitor	P154	
B Extension Lead	P157D	
C Multi input adaptor	P133G	
D Call Button	P156A	
E alternative Call Button	P156B	
F		
G		
H		
I		